



Accuray Precision® Plan Overview



Patient Name: Lin, Chih-Tung
Gender: Female
Medical ID: 1171856

Date Of Birth: 14 Jun 1976
Patient Key: 587794
Version: Accuray Precision 3.3.1.3 [2]

Plan Name: Brain SRS-7419
Plan Status: Approved
Plan Type: Standard
Treatment Machine / Serial No.: 4010594_1000MU 1180 MU/min, rev 2 / 4010594
Plan Saved Date: 14 May 2026, 1:02:21 PM (hr:min:sec) UTC+08:00

Prescription: 99.00% of PTV, 1800 cGy
Prescribed Dose per Fraction (cGy): 1800
Planned Fractions: 1
Maximum Dose (cGy): 1905
MU Planned / MU per Fraction: 14352.4 / 14352.4

Field Width (cm): 1.0, Fixed
Pitch: 0.080
Modulation Factor / Actual: 2.000 / 1.5
Sinogram Segments: 9.5
Beam On Time (sec): 825.6

Red Lasers Offset (IECf, mm): X: 0.00 Y: 20.00 Z: 0.00
Couch Travel (mm): 23.6
Couch Speed (mm/sec): 0.0
Gantry Period (sec): 29.1
Gantry Rotations: 28.4

Plan Delivery Type / Mode: Helical / IMRT
Dose Calculation Algorithm: Convolution-Superposition
Planning Method: VOLO Ultra


Dose Calculation Resolution / Type: High / Final
Spacing (IECp, mm) X: 0.88 Y: 1.00 Z: 0.88
Size (IECp, voxel) X: 604 Y: 161 Z: 528

Tracking Method: None
Image Angles (deg):
Tracking Target VOI:

Reference Point (IECp, mm): X: 30.8 Y: 5.0 Z: -10.0
Reference Dose (cGy): 1905



DVH: Absolute Dose

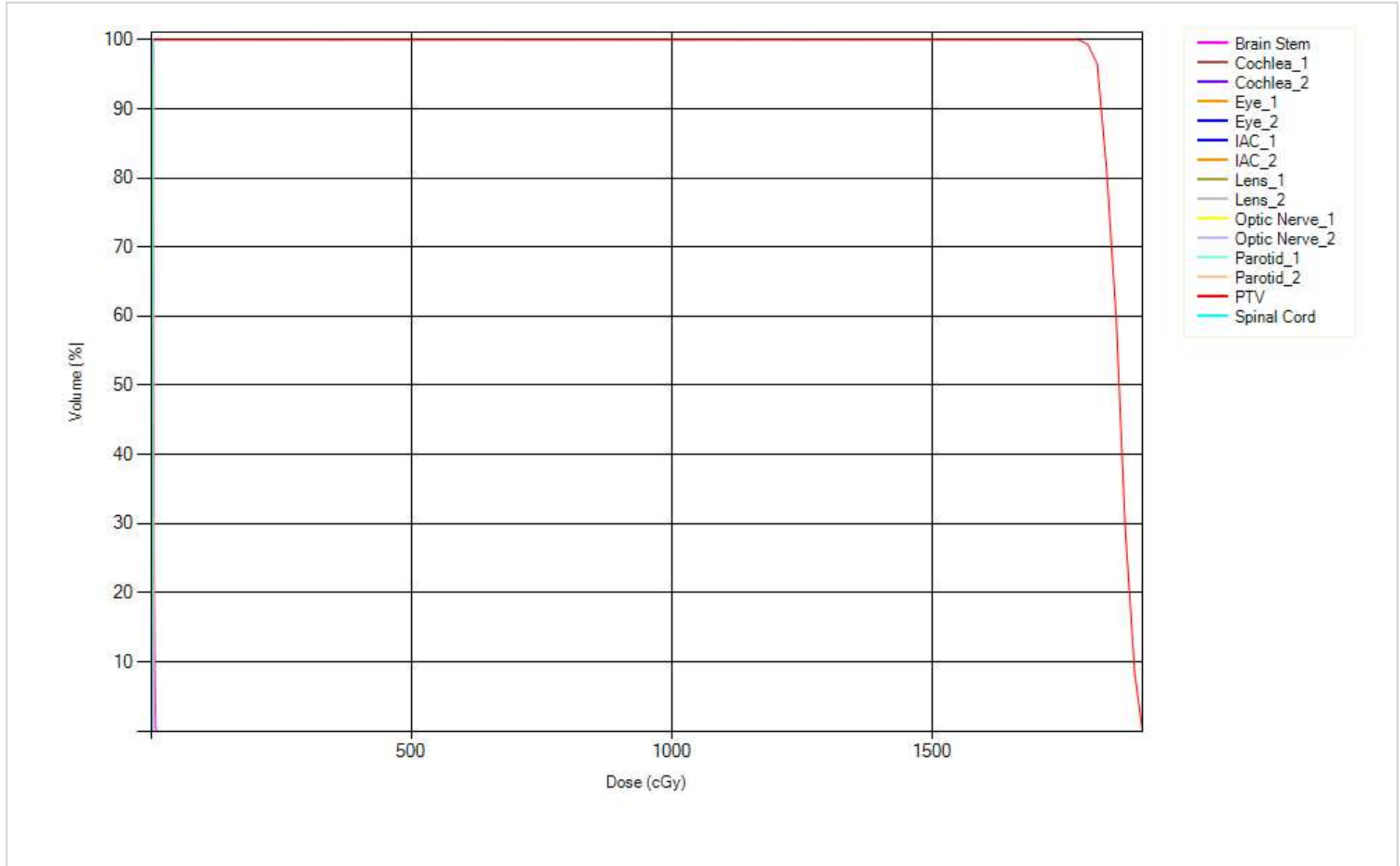



VOI

VOI List

VOI	Volume (cm³)	Min (cGy)	Mean (cGy)	Max (cGy)	CI	nCI	HI	Coverage %	Beam Inter.
Body	2559.99	0	53	1905	n/a	n/a	n/a	n/a	Allowed
Brain Stem	18.38	1	4	10	n/a	n/a	n/a	n/a	Allowed
Cochlea_1	0.40	2	3	3	n/a	n/a	n/a	n/a	Allowed
Cochlea_2	0.39	2	2	2	n/a	n/a	n/a	n/a	Allowed
Eye_1	7.99	1	2	4	n/a	n/a	n/a	n/a	Exit Only
Eye_2	8.80	1	2	4	n/a	n/a	n/a	n/a	Exit Only
Frame	3294.57	1	7	44	n/a	n/a	n/a	n/a	Allowed
IAC_1	0.38	3	3	4	n/a	n/a	n/a	n/a	Exit Only
IAC_2	0.35	2	2	3	n/a	n/a	n/a	n/a	Exit Only
Lens_1	0.34	2	2	2	n/a	n/a	n/a	n/a	Exit Only
Lens_2	0.39	2	2	2	n/a	n/a	n/a	n/a	Exit Only
Optic Nerve_1	0.65	3	4	5	n/a	n/a	n/a	n/a	Allowed
Optic Nerve_2	0.70	2	3	5	n/a	n/a	n/a	n/a	Allowed
Parotid_1	10.60	1	2	2	n/a	n/a	n/a	n/a	Allowed
Parotid_2	9.64	1	1	2	n/a	n/a	n/a	n/a	Allowed
PTV	1.78	1782	1858	1905	1.18	1.19	1.06	99.35	n/a
Spinal Cord	0.36	1	1	2	n/a	n/a	n/a	n/a	Exit Only
P01	2.93	1631	1823	1905	1.00	1.40	1.06	71.34	n/a
HOT	0.45	1812	1850	1905	4.61	4.61	1.06	100.00	n/a
RING	1.11	1782	1858	1905	1.91	1.93	1.06	98.96	n/a
All Target Regions	n/a	1631	1823	1905	0.00	0.00	0.00	0.00	n/a
All Critical Regions	n/a	1	7	44	0.00	0.00	0.00	0.00	n/a
Soft Tissue	n/a	0	52	1801	0.00	0.00	0.00	0.00	n/a



Dx Vx Values

DVH	Dose (cGy)	Dose (% of 1800 cGy Rx Dose)	Volume (cm³)	Volume (%)	Criteria	Value	
PTV	1800	100.0	1.77	99.4	Volume (%) >= 95.0, Min: 90.0	99.4	✓
Brain Stem	10	0.5	0.00	0.0	Dose (cGy) <= 600	10	✓
Cochlea_1	3	0.2	0.00	0.0	Dose (cGy) <= 400	3	✓
Eye_1	4	0.2	0.00	0.0	Dose (cGy) <= 400	4	✓
Eye_2	4	0.2	0.00	0.0	Dose (cGy) <= 400	4	✓
Lens_1	2	0.1	0.00	0.0	Dose (cGy) <= 200	2	✓
Lens_2	2	0.1	0.00	0.0	Dose (cGy) <= 200	2	✓
Optic Nerve_1	5	0.3	0.00	0.0	Dose (cGy) <= 400	5	✓
Optic Nerve_2	5	0.3	0.00	0.0	Dose (cGy) <= 400	5	✓



Density

Density Overrides

Priority	VOI Name	Relative Electron Density	Mass Density (g/cm ³)
1	Frame	n/a	0.001

Note: Density overrides are listed in priority order with the highest priority listed first. In regions where two VOIs overlap, the higher priority density override has been used.



Approval and ICD

Approved by:

Name: Lee, Yueh-Chun
User ID: Y1014
Approval Date: 14 May 2026, 1:02:34 PM (hr:min:sec) UTC08:00:00

Specialists

Radiation Oncologist: Lee, Yueh-Chun
Surgeon:
Referring Physician:
Physicist: Chang, Hank
Created in Version: Accuray Precision 3.3.1.3 [2]
Last Modified by: Chang, Hank
Saved Deliverable by: Chang, Hank
Plan Key: 589354/1
Delivery Signature: JLR7y Q55UA




Treatment Site: Secondary malignant neoplasm of brain
Treatment Site ICD: C79.31
Additional Sites:
Additional Site Codes:
Patient Position: HFS

Auto Generated Results:

Not Applicable

Name	Signature	Role	Date
		Physician	
		Physicist	

Report Snapshots

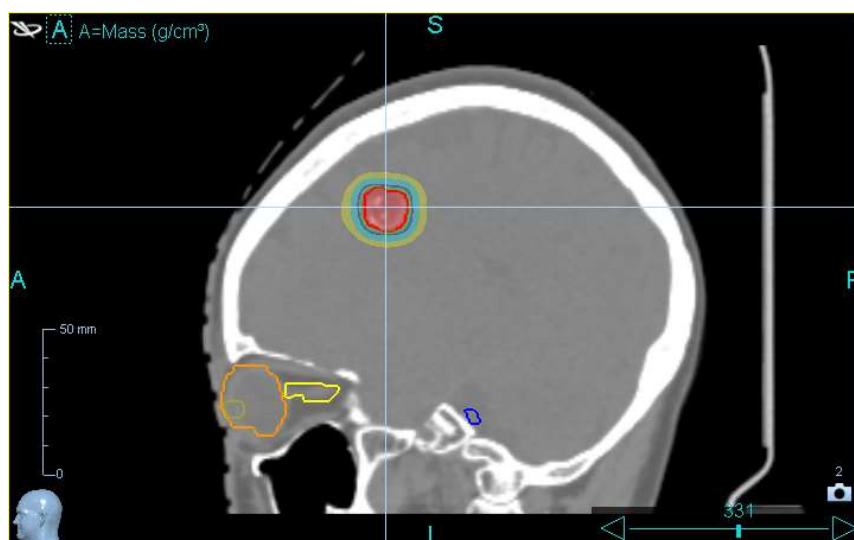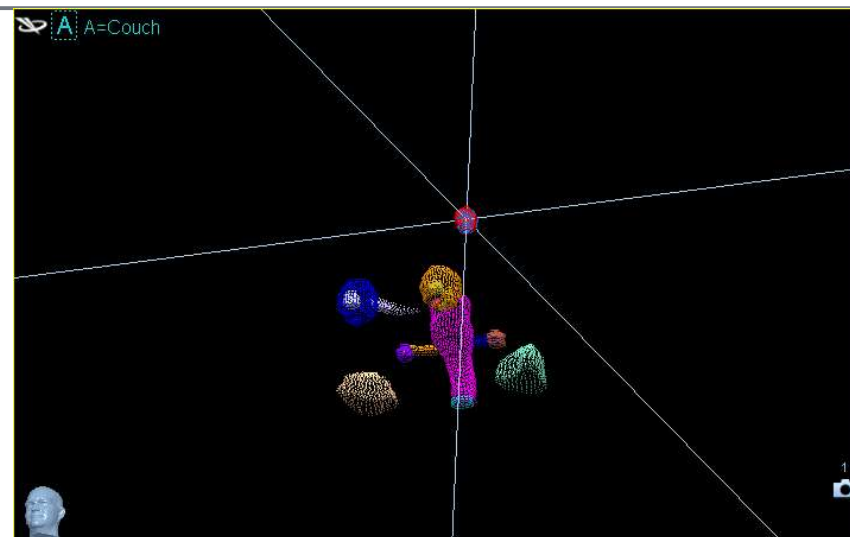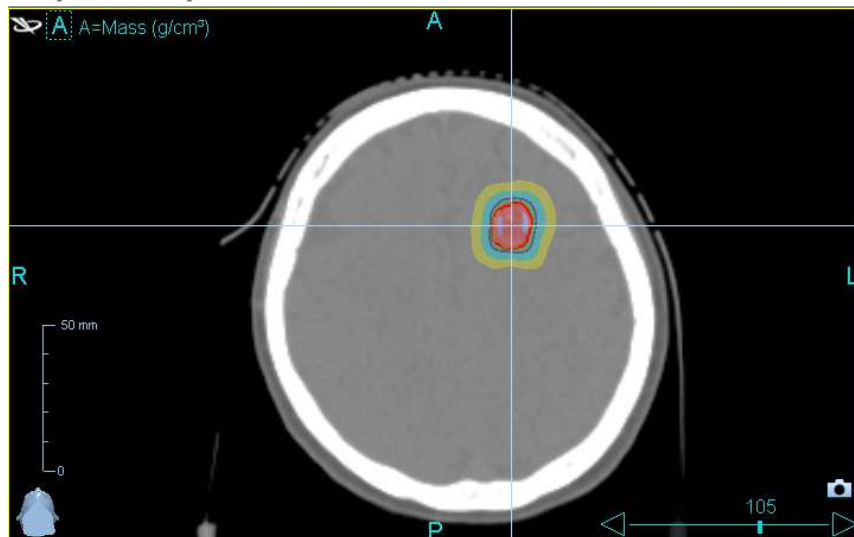



Imaging Series

Primary Series: Yes
Modality: CT
Scan Date: 14 May 2026, 09:07:30 AM (hr:min:sec)

Scanner Model: Big Bore
Patient Position: HFS
Number of Slices: 161

Series UID/Description	Study UID/Description:
1.3.46.670589.33.1.63914346626862202300002.4751642811191632899	1.3.46.670589.33.1.63914346442865800800001.5590220498408329157
	Mask SRS Brain

	X	Y	Z
Origin (IECp, mm)	-225.000	-100.000	197.400
Size (IECp, voxel)	512	161	512
Spacing (IECp, mm)	0.88	1.00	0.88